


# Energy Transfer - Wind Turbine

Model No. ET-8783

## Equipment List

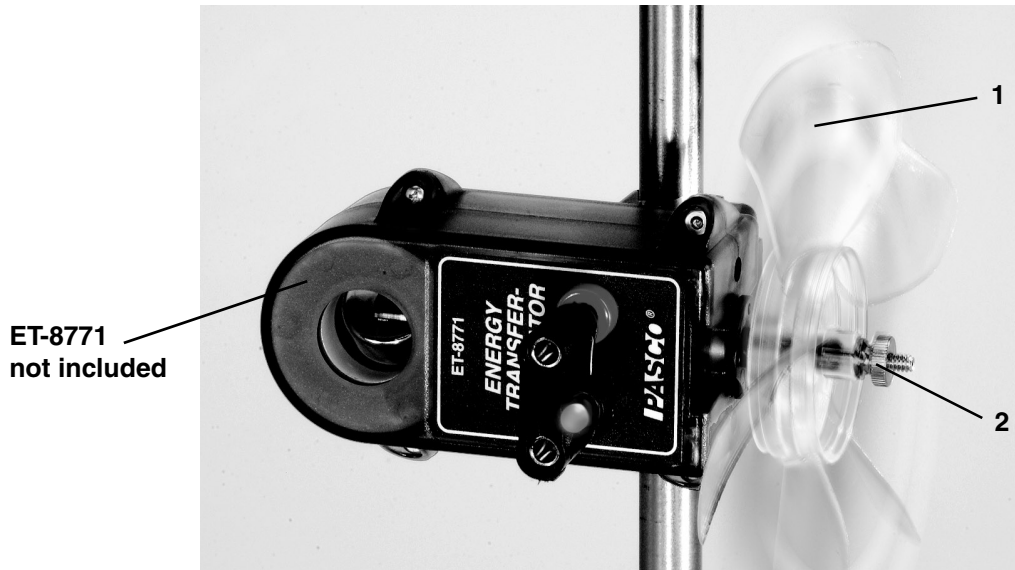

Included Equipment	Replacement Model Number*
1. 5-Blade Fan	550-023
2. Thumb Nut	614-006
3. Spacer (hidden from view)	615-046

\*Use Replacement Model Numbers to expedite replacement orders.

Additional Equipment Required	
Energy Transfer Generator	ET-8771

## Introduction

The Energy Transfer - Wind Turbine accessory can be attached to the Energy Transfer Generator (ET-8771) to create a wind energy turbine.

## Equipment Setup

### Energy Transfer-Wind Turbine Assembly

1. Remove the thumb nut and pulley from the Energy Transfer Generator. Do not discard these components. Retain for other use.
2. Place spacer over plastic shaft of the Energy Transfer-Generator.
3. Place 5-blade fan over same plastic shaft of Energy Transfer Generator.
4. Note orientation of 5-blade fan. See Figure 1.
5. Secure 5-blade fan with thumbnut. Do not over tighten.

